

Care for life Dedicate to excellence



Shanghai Kehua Bio-engineering Co., Ltd.
Mail: International@skhb.com
Phone: +86-21-64950625
Website: <https://www.skhb.com/en/>
Address: 1189 North Qinzhou Road, 200233, Shanghai

Rev.: V4 Rel.: June 6, 2024



KHB Diagnostic Kit for HIV (1+2) Antibody (Colloidal Gold) v2

WHO PQ
CE IVDR

KHB Diagnostic Kit for HIV (1+2) Antibody (Colloidal Gold) V2

Easy to Use

- Small Sample Size: 40µL
- Reads the result within 15-25 minutes
- No equipment required

Superior Raw Material

- Our HIV 1+2 rapid test who utilized superior gp 160 antigens as the major raw material, which is expressed in mammalian cells in our labs. Unlike the antigen expressed in E. Coli, this antigen has great solubility in Water, high sensitivity, and very low cross-reactivity. With this superior antigen, the sensitivity and specificity of the kit is ensured.

Stable Storage

- Storage up to 24 months at 4-30°C


Suitable for different samples

- Whole blood, serum, and plasma


Kit Components

Component (product code)	50 tests R-401-50-C-2	50 tests R-401-50-C-3	50 tests KH-R-02	50 tests KH-R-03	50 tests KH-R-02-CE	1 tests R-401-1-C-CE	25 tests R-401-25-C-1-CE	Accessories kit A-GOLD-01
HIV (1+ 2) antibody test cassette	50	50	50	50	50	1	1	N/A
Sample diluent	4ml x 1 bottle	4ml x 2 bottle	4ml x 1 bottle	4ml x 2 bottle	4ml x 1 bottle	0.5ml x 1 bottle	4ml x 1 bottle	N/A
Safety lancet	50	50	N/A	N/A	N/A	1	25	50
Alcohol pad	50	50	N/A	N/A	N/A	1	25	50
Disposable transfer pipette	50	50	N/A	N/A	N/A	1	25	50
Instruction for use	1 pcs	1 pcs	1 pcs	1 pcs	1 pcs	1 pcs	1 pcs	1 pcs

Test Procedure

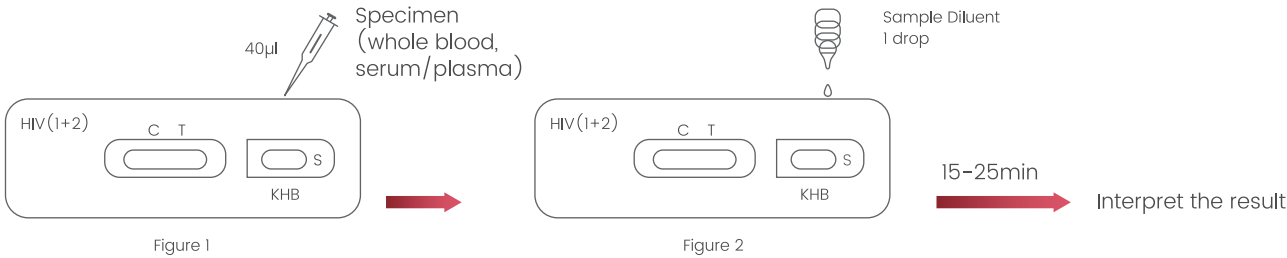

Result Analysis

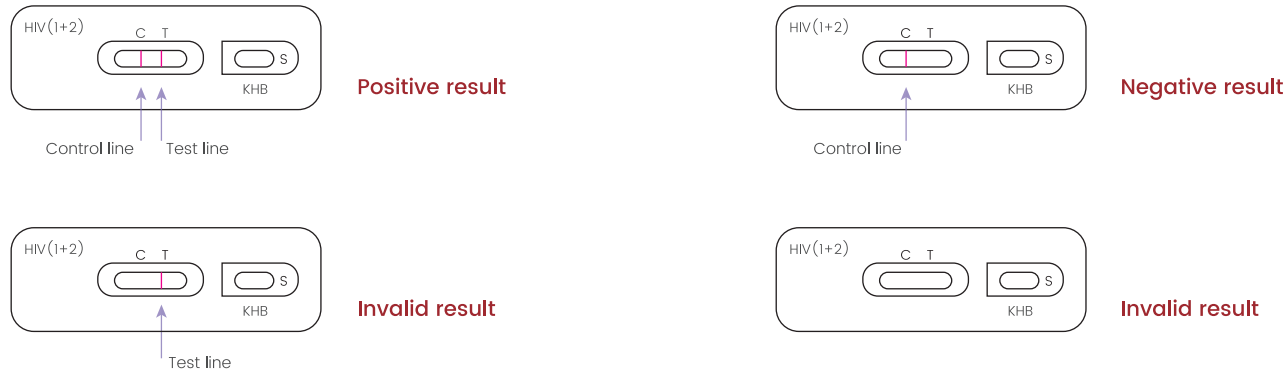

Performance features

	Initial(95% CI)	Final(95% CI)
Sensitivity%	100% (99.2%-100%)	100% (99.2%-100%)
Specificity%	100% (99.4% -100%)	100% (99.4% -100%)
Invalid rate%	0%	
Inter- reader variability%	0.09%	

	Reference testing results			
Diagnostic Kit for HIV (1+2) Antibody (Colloida Gold) v2 Final results	HIV +	HIV -	Total	
	HIV (1&2) +	460	0	460
	HIV-1+	0	0	0
	HIV-2+	0	0	0
	HIV-	0	657	657
	Total	460	657	1117

- Sensitivity (95% CI): 460/460=100% (99.2% -100%)
- Specificity (95% CI): 657/657=100% (99.4%-100%)
- Indeterminate results: 0/1117=0%
- False negative results: 0/657=0% (See IFU Annex3)
- False positive results: 0/460=0% (See IFU Annex3)